



Barrier Technologies®

Your Safety is Our Concern™

Radiation Protection Accessories



Mobile Barriers

Barrier Technologies® Mobile Barriers are made with **Barrier Safe™** Radiation Protection Glass — a laminated, safety, radiation shielding lead glass with **2.0mm Pb equivalent @ 150 kVp** to provide maximum protection.

- Safety - Shatter Resistant Glass Exceeds ANSI Z97.1-1984 & 16 CFR 1201 Cat.II requirements.
- Will not discolor
- Delivered fully assembled



Advantages Over Lead Acrylic:

- Much lighter & easier to maneuver
- Better light transmission
- Superior visual clarity
- Greater resistance to scratching

Corning® Med-X® Glass

APRON RACKS



Spider Rack

- Porta-Rack stable apron rack
- Hangers swing forward and backwards
- Ideal for OR and Cath Labs
- 24" Wheel Base
- Available with top mounted glove rack



Wall Mount Swing Arm Apron Rack

- Lead apron rack
- Space saving wall console
- 5 hangers
- Available in right or left swing



Mobile Swing Arm Rack

- Stable mobile apron rack
- Solid steel construction
- 25"W x 28"L x 55"H



Deluxe Hanger



Peg Rack

- Strong steel construction
- White anodized finish
- Easy to install
- Available with top mounted glove rack


Valet Rack

- Heavy duty steel construction
- Holds up to 10 aprons
- 4 Screws for easy install



Table Shielding Drapes

Custom Table Scatter Shield features a saddle back design allowing a secure fit to any style imaging table. The overlapping paneled shields are secured to a padded surface that simply drapes over the table-top, beneath the mattress.

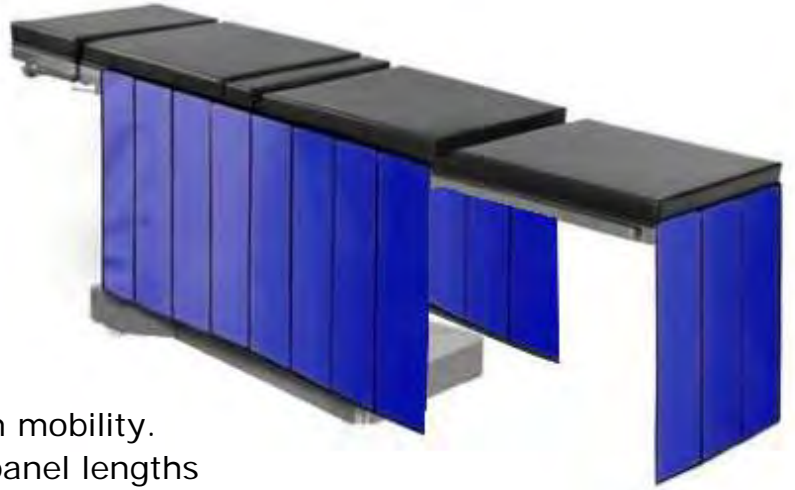

- Pleated lead panels enable staff and c-arm mobility.
- Customized to specific table dimensions, panel lengths or to include table end flaps.
- Constructed of 0.5mm Pb equivalent panels.
- Select from our wide range of UltraFlex[™] anti-bacterial/stain resistant fabrics

Protective Blankets

Nylon-covered 0.5mm Pb equivalent protective blanket may be used for shielding patients or equipment.



Porta Shield

Porta-Shield is ideal for patient protection during diagnostic procedures and features adjustable height with sturdy five leg base for stable mobility. Available in three sizes and covered in Royal Blue UltraFlex[™].

INTRODUCING



X-RAY APRON CLEANING WIPES

COVID-19
COMPLIANT


 **Barrier Technologies**[®]
Your Safety is Our Concern[™]



**COMPREHENSIVELY ACTIVE
AGAINST STAINS AND BACTERIA**



WIPE DOWN APRON

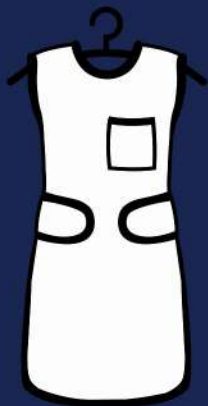

HANG AND AIR DRY



KILL HARMFUL BACTERIA

www.BarrierTechnologies.com



TERMINATOR XR™

MAXIMUM RADIATION PROTECTION

1.0mm Pb LE Protection Per Panel
2.0mm Pb LE Protection In Overlapping Areas
100% attenuation at 160 kVp




ULTRAFLEX™

Covered in UltraFlex™
Anti-bacterial stain resistant fabric.

- Body Contoured
- Height Adjustable
- Completely Mobile
- Fully Encapsulating



Enhanced Operator Comfort & Safety

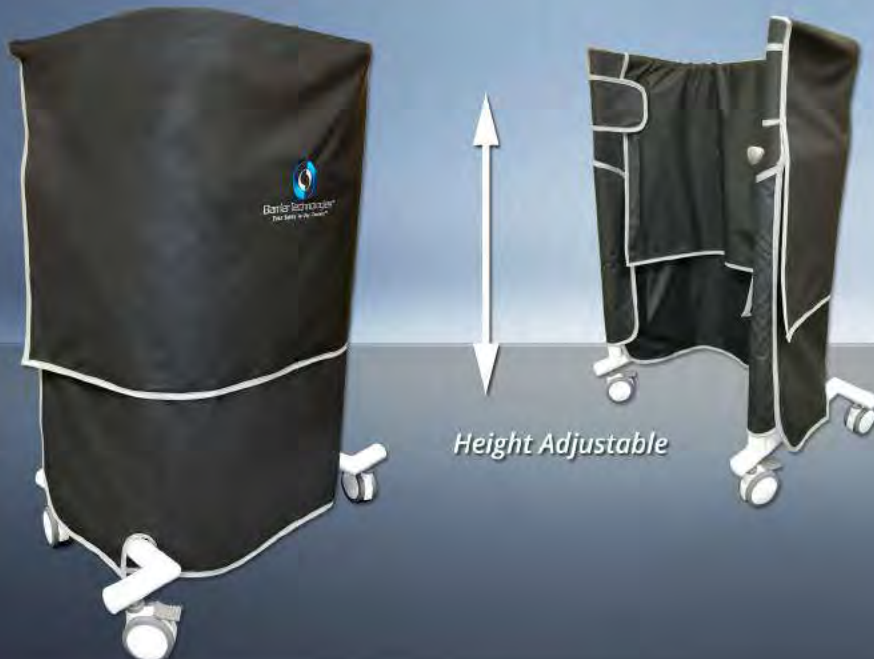

Height Adjustable

- Virtual Weightlessness
- Unimpeded Patient Access
- Optimal Scatter Protection
- Reduces Radiation Exposure

See All of Our Radiation Protection Products at:
WWW.BARRIERTECHNOLOGIES.COM



Aprons



Leaded Eyewear



Gloves



Thyroid Collars



Mobile Barriers



Accessories



Secure Shields



Dental



7060 W. STATE ROAD 84, Suite 8, DAVIE, FL 33317

Tel: (800) 353-4350 ▪ Fax: (800) 353-4351

Tel: (954) 921-0330 ▪ Fax: (954) 921-0331

WWW.BARRIERTECHNOLOGIES.COM